

Thermo Scientific Finnpiquette C1

The Finnpiquette C1

This portable, lightweight serological pipette holder is designed for use with graduated and volumetric glass and plastic pipets from 1 ml to 100 ml. The Finnpiquette C1 offers simple, efficient pipetting performance with maximum pipetting comfort.

Total control of pipetting


The Finnpiquette C1 excels in offering total control of pipetting due to eight speed alternatives for both aspiration and dispensing. It also enables two different ways of “gravity only” dispensing. Additional control is gained from the 1 ml pipette support that stabilizes pipetting with narrow 1 ml pipets. The Finnpiquette C1 also has a large illuminated LC display, that shows battery charge and pump speeds for both aspiration and dispensing. The display alarms with a blinking light when the battery charge gets low.

Lithium-ion battery technology

Effective and comfortable to operate, the Finnpiquette C1 features a heavy-duty lithium-ion battery that allows the pipettor to operate continuously, approximately three times longer than regular units using nickel cadmium batteries. The lithium-ion battery can be recharged in just three hours. A universal power supply is automatically delivered with the unit.

Compact storage between applications

Both a wall-mounted holder and a table-mounted holder are included in the package.

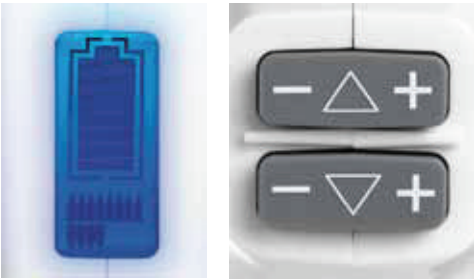

- Lithium-ion technology – provides up to 15 hours of continuous operating
- Lightweight ergonomic design
- Autoclavable silicon adapter
- 1 ml pipette support
- Two ways of “gravity only” dispensing
- Universal power supply
- Fast recharging
- 2-year warranty with web registration

Finnpiquette C1 Technical Specifications

Aspirate Speeds	8 options
Dispense Speeds	8 options (+ gravity dispense)
Battery	Lithium-ion
Approx. Runtime Between Charges	15 hours of continuous use
Recharge Time	3 hours
Weight	220 g
Nosecone Pipet Holder	Autoclavable silicone rubber
Pipet Compatibility	All major brands of glass or plastic serological pipets, 1 ml to 100 ml
Filter Assembly	0.45 micron standard
Filter Material	PTFE

Product specifications are subject to change without prior notice.

Ordering Information

Cat. no.	Description	Qty / Unit
4580800	Finnpiquette C1	1
Filters		
4580560	Filter, non-sterile	5/bag
4580570	Filter, sterile	5/bag
Accessories and spare parts		
22280	Pipette gripper housing	1
22278	Silicone pipette gripper	1
22300	Table mounted holder	1
22302	Wall mounted holder	1